

reaClean professional

Our reaClean professional Wash Check indicators were developed to meet the ever-increasing requirements of hygiene monitoring in medicine.

Various studies show that reaching the specified temperatures is decisive for the effective elimination of germs such as viruses, bacteria and fungi. If, for example, the surface temperature remains below 85 °C (min. 16 min) during a cleaning process, numerous germs including the persistent norovirus can survive.



The two wash indicators (active and passive) document the cleaning parameters such as time, strength of the cleaning solution and mechanical action of the water during the cleaning process. To control and document the cleaning and disinfection of medical instruments even more reliably, we have added our reaPoint professional (88°C or 93°C) for precise temperature control.

Our reaPoint professional are irreversible temperature indicators predestined for the control and documentation of the temperature during the cleaning process. It indicates by a color change from white to black that the specified temperature has been reached.

The additional temperature monitoring is done cost-effectively by a single reaClean professional indicator (complies with the specifications from ISO 15883-5). The indicator will be placed in a supplied holder in the machine.

REATEC AG


Phone: +41 525 880 298 Mail: office@reatec.ch

Klettgauerstr. 67 CH-8212 Neuhausen am Rheinfall

reaClean professional



product example: 088.CPX.03.XXX

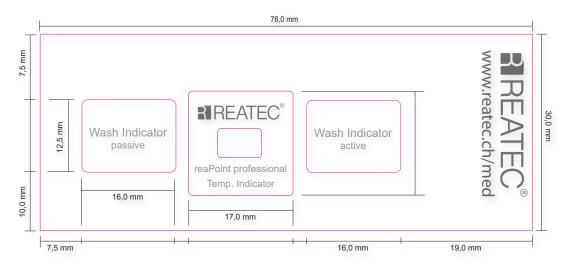

product number	Temp. °C	
088.CPX.03.001.XXX	88	Starter Kit professional 88°C / 250 pcs
093.CPX.03.001.XXX	93	Starter Kit professional 93°C / 250 pcs
XXX.CSX.02.002.XXX ¹	-	Starter Kit standard / 250 pcs
088.CPX.03.001.0250	88	Refill pack professional 88°C / 250 pcs
093.CPX.03.001.0250	93	Refill pack professional 93°C / 250 pcs
XXX.CSX.02.002.0250 ¹	-	Refill pack standard / 250 pcs
XXX.CHF.XX.003.0001	-	Wash Check holder

¹ standard wash check without temperature measuring point.

The Starter Kit	comes in a h	nigh auglity	hav and	contains.

- 1 reaClean stainless steel holder
- 250 reaClean wash check indicators
- 250 reaPoint professional temperature measuring points (only for Starter Kit professional)
- instruction manual for use of the indicator
- documentation manual for recording the determined results

Tyvek paper, PET (plastic)	material
2-3	measuring points
yes	protected
yes	dishwasher suitable
18 months after opening	shelf life
whole wash cycle	reaction time wash check
1-2 sec.	reaction time temperature
76 x 30 mm	dimensions
cool, dry	storage conditions